


436 / 540 nm

Method(e) / Método


0721

pH 6.5 - 8.2


### Nullwert / Blanc value / Zéro / Cero

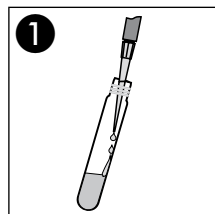


10 ml Probe  
Sample  
Enchantillon  
Muestra

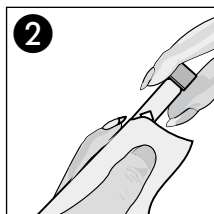

Säubern  
Clean  
Nettoyer  
Limpiar

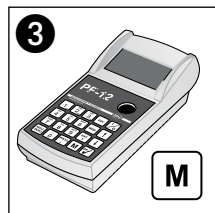


Messung  
Measurement  
Mesure  
Medición

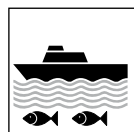


Meerwasser  
Sea water  
Eau de mer  
Agua de mar


### Messwert / Sample / Echantillon / Muestra

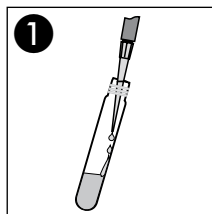


10 ml Probe  
Sample  
Enchantillon  
Muestra

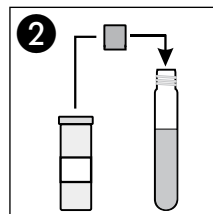


1 x NANOFIX pH

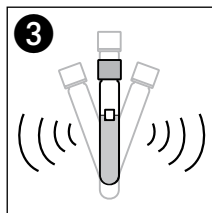


Schütteln  
Shake  
Agiter  
Agitar

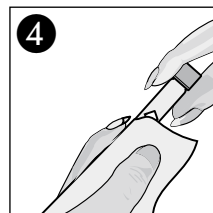

Säubern  
Clean  
Nettoyer  
Limpiar

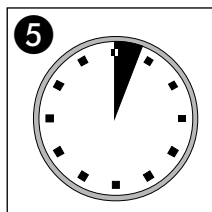

3'00 min

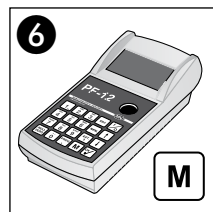

Messung  
Measurement  
Mesure  
Medición