

NANOCOLOR® Ammonium 10

Amonio


690 nm

Method(e) / Método


0041	0.2 - 8.0 mg/l $\text{NH}_4\text{-N}$
0042	0.2 - 10.0 mg/l NH_4^+
0043	0.2 - 10.0 mg/l NH_3

Test 0-04


REF 985 004

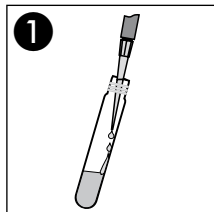


1
1.0 ml Probe
Sample
Enchantillon
Muestra

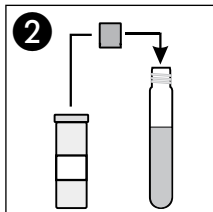


2
1 x **NANOFIX R2**

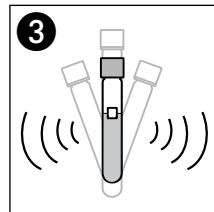


3
Schütteln
Shake
Agiter
Agitar

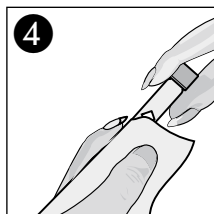


4
Säubern
Clean
Nettoyer
Limpiar

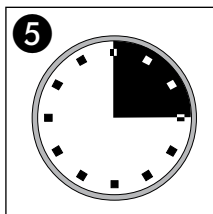

5
15'00 min

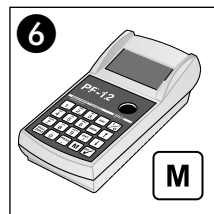

6
Messung
Measurement
Mesure
Medición

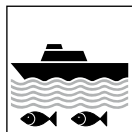

Meerwasser
Sea water
Eau de mer
Agua de mar