

NANOCOLOR® Ammonium 200

Amonio


585 nm

Method(e) / Método


0061	30 - 160 mg/l $\text{NH}_4\text{-N}$
0062	40 - 200 mg/l NH_4^+
0063	40 - 200 mg/l NH_3

Test 0-06


REF 985 006

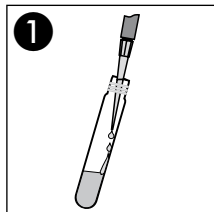


200 µl Probe
Sample
Enchantillon
Muestra

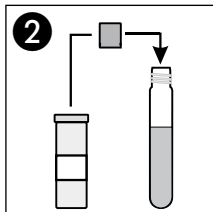


1 x NANOFIX R2

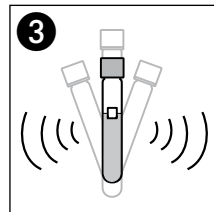


Schütteln
Shake
Agiter
Agitar

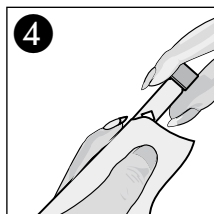


Säubern
Clean
Nettoyer
Limpiar

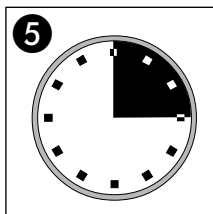

15'00 min

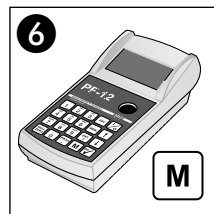

Messung
Measurement
Mesure
Medición

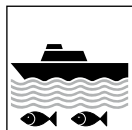

Meerwasser
Sea water
Eau de mer
Agua de mar