

NANOCOLOR® Ammonium 100

Amonio


585 nm

Method(e) / Método


0081	4 - 80 mg/l $\text{NH}_4\text{-N}$
0082	5 - 100 mg/l NH_4^+
0083	5 - 100 mg/l NH_3

Test 0-08


REF 985 008

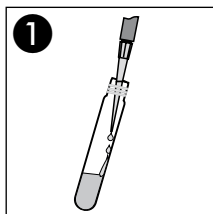


1
500 µl Probe
Sample
Enchantillon
Muestra

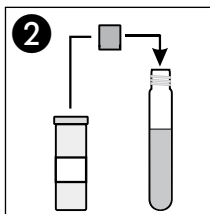


2
1 x NANOFIX R2

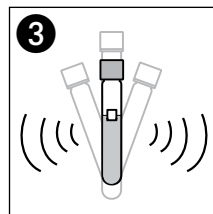


3
Schütteln
Shake
Agiter
Agitar

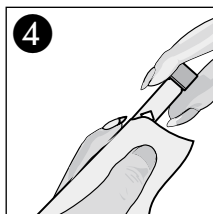


4
Säubern
Clean
Nettoyer
Limpiar

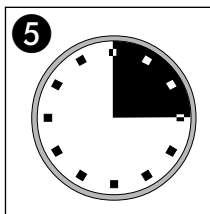

5
15'00 min

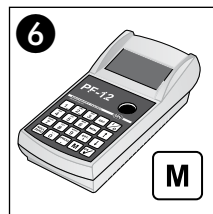

6
Messung
Measurement
Mesure
Medición

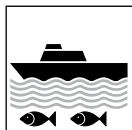

Meerwasser
Sea water
Eau de mer
Agua de mar