

# NANOCOLOR® CSB 60 000

COD / DCO / DQO

620 nm


Method(e) / Método

0121


5.0 - 60.0 g/l O<sub>2</sub> (5000 - 60000 mg/l O<sub>2</sub>)

**Test 0-12**


REF 985 012

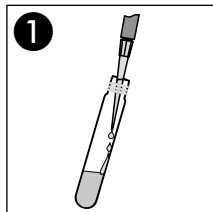


**200 µl** Probe  
Sample  
Enchantillon  
Muestra

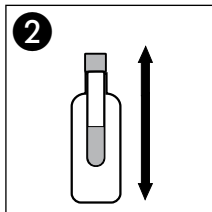


Sicherheitsgefäß  
Safety bottle  
Récipient de sécurité  
Recipiente de seguridad

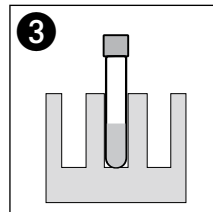


**148 °C / 2 h**

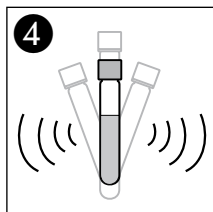


Schütteln  
Shake  
Agiter  
Agitar

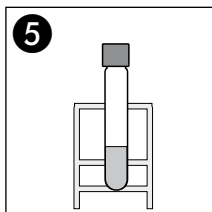

Abkühlen lassen  
Allow to cool  
Laisser refroidir  
Dejar enfriar

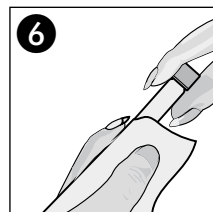

Säubern  
Clean  
Nettoyer  
Limpiar

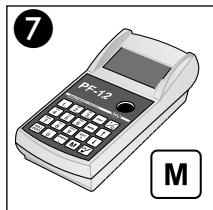

Messung  
Measurement  
Mesure  
Medición

Meerwasser  
Sea water  
Eau de mer  
Agua de mar

