

NANOCOLOR® Cadmium 2

Cadmio

Test 0-14

REF 985 014


540 nm

Method(e) / Método


0141

0.05 - 2.00 mg/l Cd²⁺


Nullwert / Blanc value / Zéro / Cero

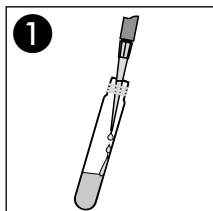


4.0 ml dest. Wasser
dist. water
d'eau distillée
agua dest.

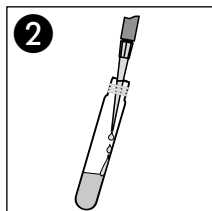


200 µl R2

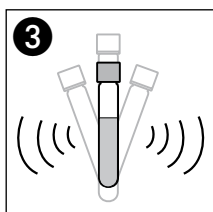


Schütteln
Shake
Agiter
Agitar

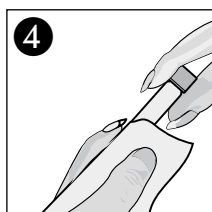

Säubern
Clean
Nettoyer
Limpiar

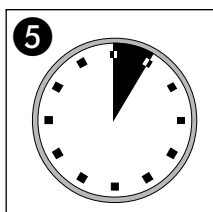

5'00 min

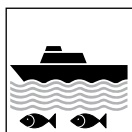


Meerwasser / Sea water
Eau de mer / Agua de mar




Messung
Measurement
Mesure
Medición


Messwert / Sample / Echantillon / Muestra

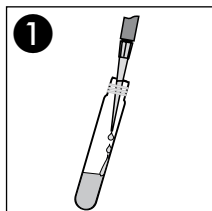


4.0 ml Probe
Sample
Echantillon
Muestra

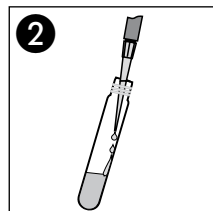


200 µl R2

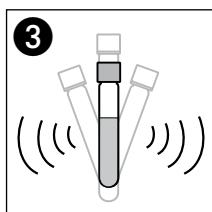


Schütteln
Shake
Agiter
Agitar

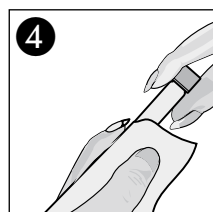

Säubern
Clean
Nettoyer
Limpiar

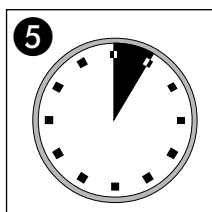

5'00 min

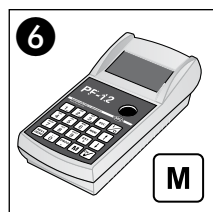

Messung
Measurement
Mesure
Medición