

# NANOCOLOR® Chlordioxid 5

Chlorine dioxide / Bioxyde de chlore /  
Dioxido de cloro

540 nm


Method(e) / Método

0181


0.15 - 5.00 mg/l ClO<sub>2</sub>

**Test 0-18**


REF 985 018

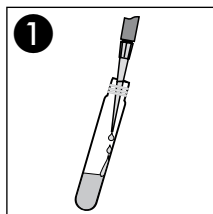


**1**  
4.0 ml Probe  
Sample  
Enchantillon  
Muestra

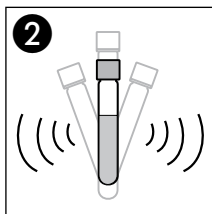


**2**  
Schütteln  
Shake  
Agitear  
Agitar

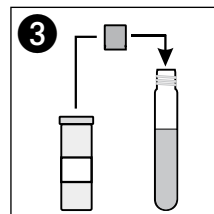


**3**  
1 x NANOFIX R2

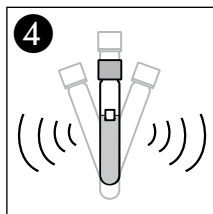


**4**  
Schütteln  
Shake  
Agitear  
Agitar

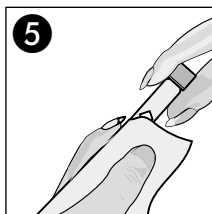

**5**  
Säubern  
Clean  
Nettoyer  
Limpiar

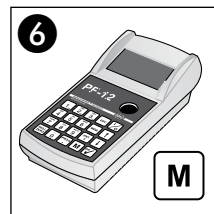

**6**  
Messung  
Measurement  
Mesure  
Medición

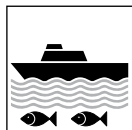

Meerwasser  
Sea water  
Eau de mer  
Agua de mar