

NANOCOLOR® Chlorid 200

Chloride / Chlorure / Cloruro

470 nm


Method(e) / Método

0191


5 - 200 mg/l Cl⁻

Test 0-19


REF 985 019

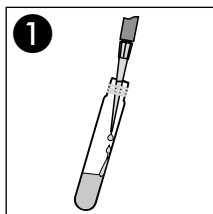


1
1.0 ml Probe
Sample
Enchantillon
Muestra

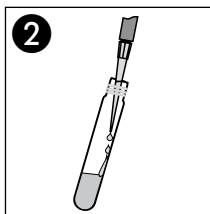


2
1.0 ml R2

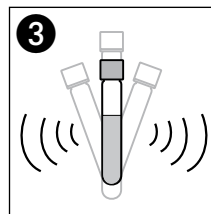


3
Schütteln
Shake
Agiter
Agitar

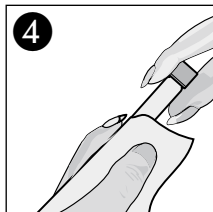


4
Säubern
Clean
Nettoyer
Limpiar

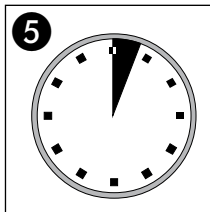

5
3'00 min

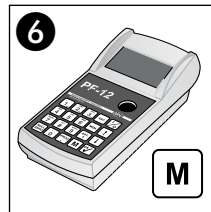

6
Messung
Measurement
Mesure
Medición

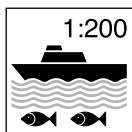

1:200
Meerwasser
Sea water
Eau de mer
Agua de mar