

NANOCOLOR® Chlorid 50

Chloride / Chlorure / Cloruro

470 nm


Method(e) / Método

0211


0.5 - 50.0 mg/l Cl⁻

Test 0-21


REF 985 021

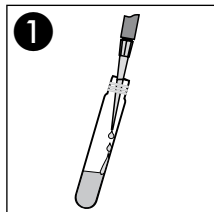


1
4.0 ml Probe
Sample
Enchantillon
Muestra

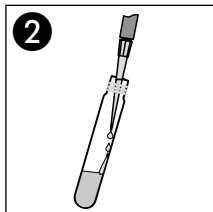


2
1.0 ml R2

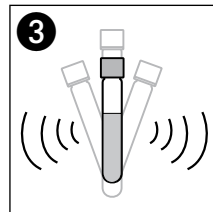


3
Schütteln
Shake
Agiter
Agitar

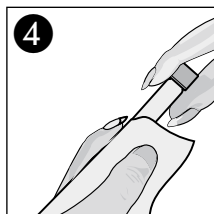

4
Säubern
Clean
Nettoyer
Limpiar

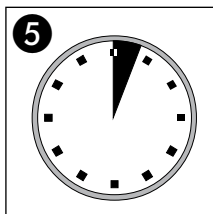

5
3'00 min

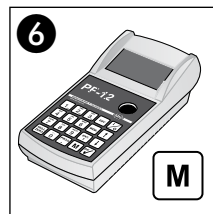

6
Messung
Measurement
Mesure
Medición



Meerwasser
Sea water
Eau de mer
Agua de mar