

NANOCOLOR® Chromat 5

Chromate / Cromat

Test 0-24


REF 985 024

540 nm


Method(e) / Método

0241 0.05 - 2.00 mg/l Cr(VI)


0242 0.1 - 4.0 mg/l CrO₄²⁻

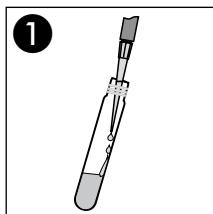


4.0 ml Probe
Sample
Enchantillon
Muestra

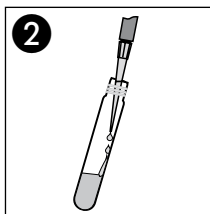


200 µl R2

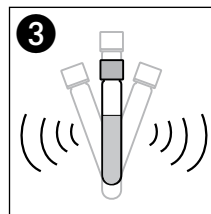


Schütteln
Shake
Agiter
Agitar

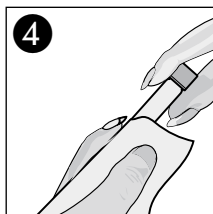


Säubern
Clean
Nettoyer
Limpiar

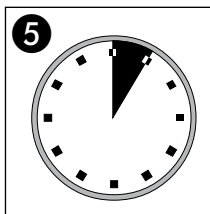

5'00 min

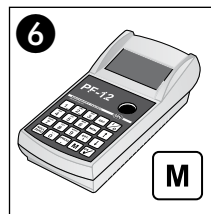

Messung
Measurement
Mesure
Medición

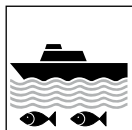

Meerwasser
Sea water
Eau de mer
Agua de mar