

NANOCOLOR® CSB 40

COD / DCO / DQO

345 nm


Method(e) / Método

0271


2 - 40 mg/l O₂

Test 0-27


REF 985 027

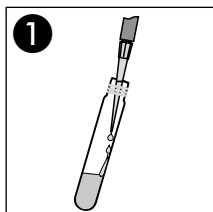


2.0 ml Probe
Sample
Enchantillon
Muestra

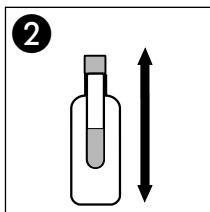


Sicherheitsgefäß
Safety bottle
Récipient de sécurité
Recipiente de seguridad

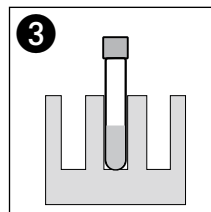


148 °C / 2 h

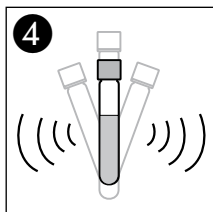


Schütteln
Shake
Agiter
Agitar

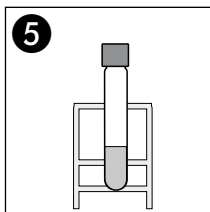

Abkühlen lassen
Allow to cool
Laisser refroidir
Dejar enfriar

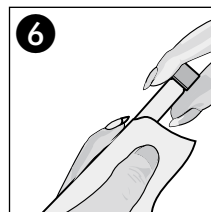

Säubern
Clean
Nettoyer
Limpiar

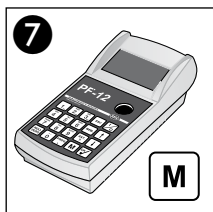

Messung
Measurement
Measure
Medición

Meerwasser
Sea water
Eau de mer
Agua de mar

