

# NANOCOLOR® CSB 15 000

COD / DCO / DQO

620 nm


Method(e) / Método

0281


1.0 - 15.0 g/l O<sub>2</sub> (1000 - 15000 mg/l O<sub>2</sub>)

**Test 0-28**


REF 985 028

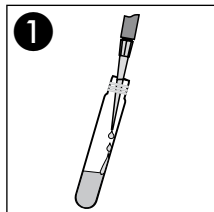


**1**  
200 µl Probe  
Sample  
Enchantillon  
Muestra

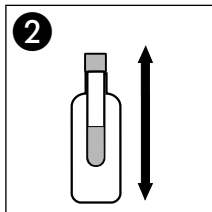


**2**  
Sicherheitsgefäß  
Safety bottle  
Récipient de sécurité  
Recipiente de seguridad

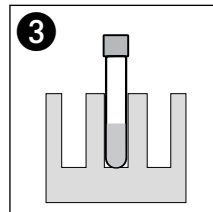


**3**  
148 °C / 2 h

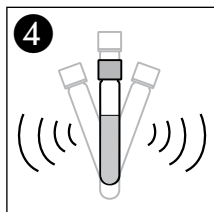


**4**  
Schütteln  
Shake  
Agiter  
Agitar

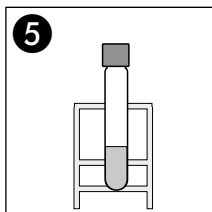

**5**  
Abkühlen lassen  
Allow to cool  
Laisser refroidir  
Dejar enfriar

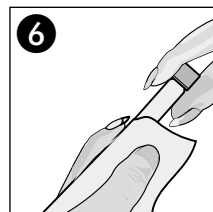

**6**  
Säubern  
Clean  
Nettoyer  
Limpiar

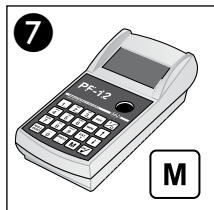

**7**  
Messung  
Measurement  
Mesure  
Medición

Meerwasser  
Sea water  
Eau de mer  
Agua de mar

