

NANOCOLOR® CSB 1500

COD / DCO / DQO

620 nm


Method(e) / Método

0291


100 - 1500 mg/l O₂

Test 0-29


REF 985 029

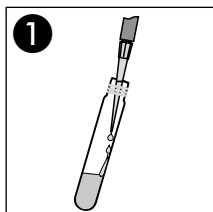


1
2.0 ml Probe
Sample
Enchantillon
Muestra

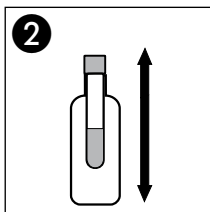


2
Sicherheitsgefäß
Safety bottle
Récipient de sécurité
Recipiente de seguridad

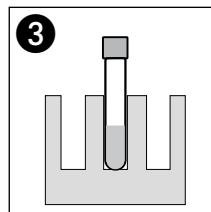


3
148 °C / 2 h

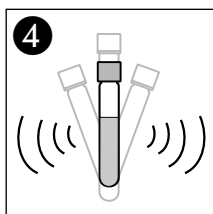


4
Schütteln
Shake
Agiter
Agitar

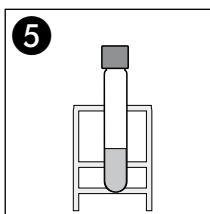

5
Abkühlen lassen
Allow to cool
Laisser refroidir
Dejar enfriar

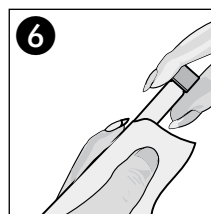

6
Säubern
Clean
Nettoyer
Limpiar

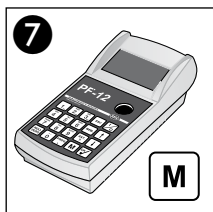

7
Messung
Measurement
Measure
Medición

Meerwasser
Sea water
Eau de mer
Agua de mar

