

# NANOCOLOR® Cyanid 08

Cyanid / Cyanure / Cianuro

585 nm


Method(e) / Método

0311


0.02 - 0.80 mg/l CN<sup>-</sup>

**Test 0-31**


REF 985 031

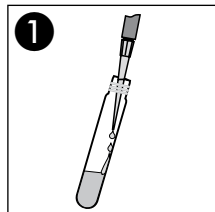


**1**  
4.0 ml Probe  
Sample  
Enchantillon  
Muestra

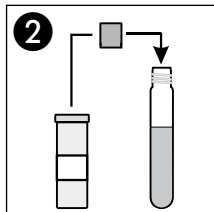


**2**  
1 x **NANOFIX R2**

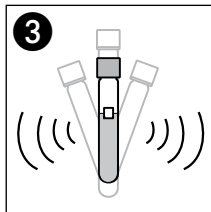


**3**  
Schütteln  
Shake  
Agiter  
Agitar

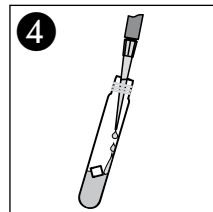


**4**  
500 µl R3

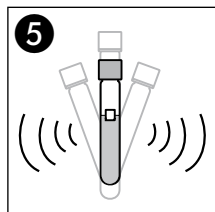


**5**  
Schütteln  
Shake  
Agiter  
Agitar

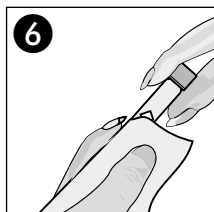


**6**  
Säubern  
Clean  
Nettoyer  
Limpiar

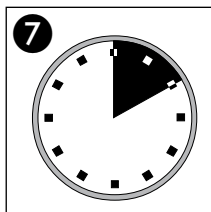

**7**  
10'00 min

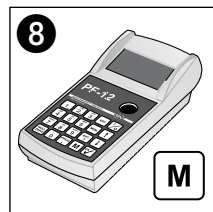

**8**  
Messung  
Measurement  
Mesure  
Medición

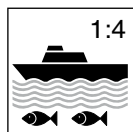

1:4  
Meerwasser  
Sea water  
Eau de mer  
Agua de mar