

NANOCOLOR® Anionische Tenside 4

Test 0-32

Anionic surfactants / Tensio-actifs
anioniques / Tensioactivos aniónicos


REF 985 032

620 nm


Method(e) / Método

0321


0.20 - 4.00 mg/l MBAS

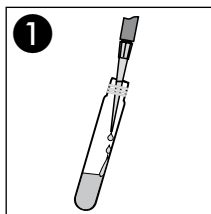


4.0 ml Probe
Sample
Enchantillon
Muestra

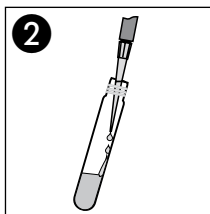


500 µl R2

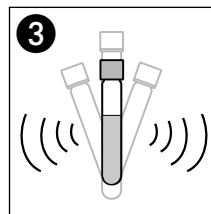


1 min Schütteln
Shake
Agiter
Agitar

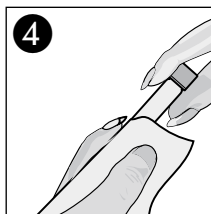


Säubern
Clean
Nettoyer
Limpiar

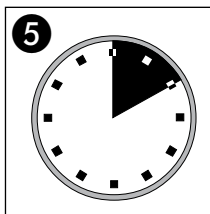

10'00 min
Phasentrennung abwarten
wait for phase separation
laisser séparer les phases
esperar la separación de fases

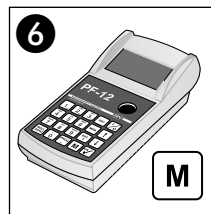

Messung
Measurement
Mesure
Medición

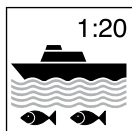

1:20
Meerwasser
Sea water
Eau de mer
Agua de mar