

NANOCOLOR® DEHA 1

Test 0-35


REF 985 035

540 nm


Method(e) / Método

0351


0.05 - 1.00 mg/l DEHA

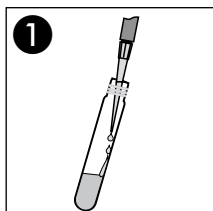


4.0 ml Probe
Sample
Enchantillon
Muestra

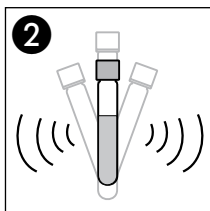


Schütteln
Shake
Agiter
Agitar

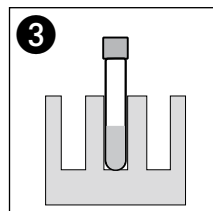


100 °C / 15 min

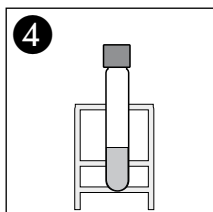


Abkühlen lassen
Allow to cool
Laisser refroidir
Dejar enfriar

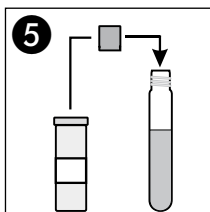


1 x NANOFIX R2

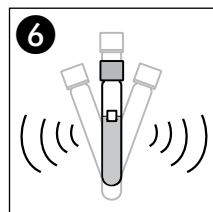


Schütteln
Shake
Agiter
Agitar

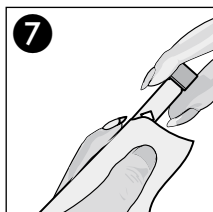


Säubern
Clean
Nettoyer
Limpiar

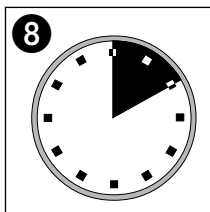

10'00 min

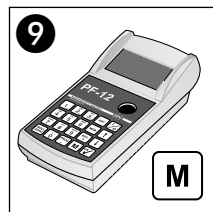

Messung
Measurement
Mesure
Medición

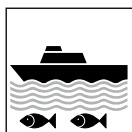

Meerwasser
Sea water
Eau de mer
Agua de mar