

NANOCOLOR® Eisen 3

Iron / Fer / Hierro

540 nm


Method(e) / Método

0371


0.10 - 3.00 mg/l Fe

Test 0-37


REF 985 037

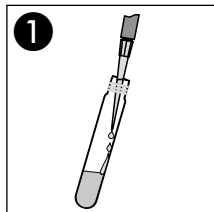


1
4.0 ml Probe
Sample
Enchantillon
Muestra

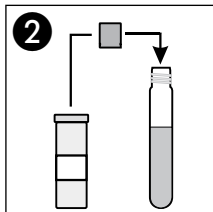


2
1 x **NANOFIX R2**

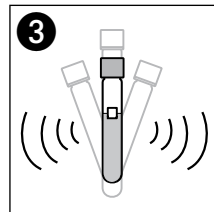


3
Schütteln
Shake
Agiter
Agitar

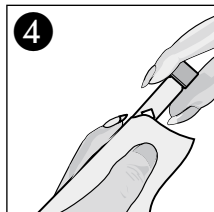


4
Säubern
Clean
Nettoyer
Limpiar

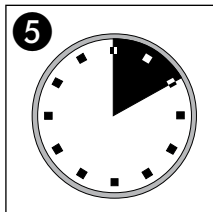

5
10'00 min

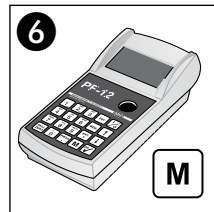

6
Messung
Measurement
Mesure
Medición

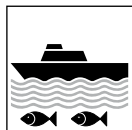

Meerwasser
Sea water
Eau de mer
Agua de mar