

# NANOCOLOR® Fluorid 2

Fluoride / Fluorure / Fluoruro

620 nm


Method(e) / Método

0401


0.1 - 2.0 mg/l F<sup>-</sup>

**Test 0-40**


REF 985 040

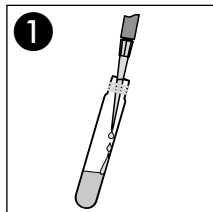


**1**  
2.0 ml Probe  
Sample  
Enchantillon  
Muestra

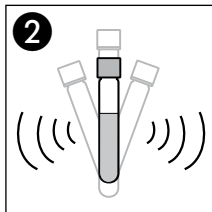


**2**  
Schütteln  
Shake  
Agiter  
Agitar

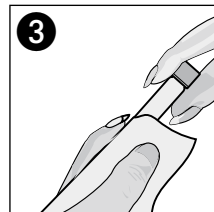


**3**  
Säubern  
Clean  
Nettoyer  
Limpiar

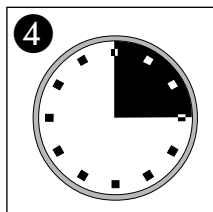

**4**  
15'00 min

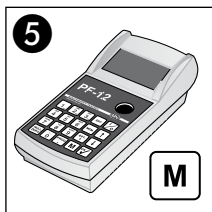

**5**  
Messung  
Measurement  
Mesure  
Medición

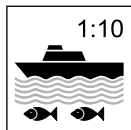

1:10  
Meerwasser / Sea water  
Eau de mer / Agua de mar