

NANOCOLOR® Kalium 50

Potassium / Potasio

690 nm


Method(e) / Método

0451


2 - 50 mg/l K⁺

Test 0-45


REF 985 045

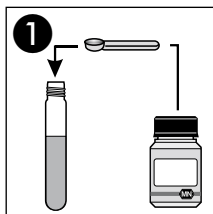


1 x R2

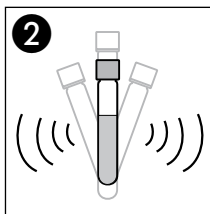


Schütteln
Shake
Agiter
Agitar

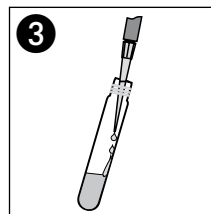


2.0 ml Probe
Sample
Enchantillon
Muestra

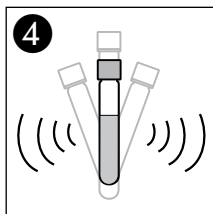


10 s Schütteln
Shake
Agiter
Agitar

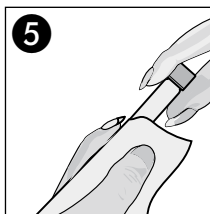

Säubern
Clean
Nettoyer
Limpiar

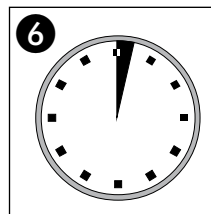

2'00 min

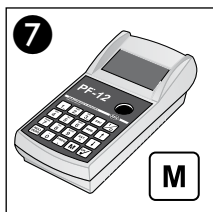

Messung
Measurement
Measure
Medición

Meerwasser
Sea water
Eau de mer
Agua de mar

