

NANOCOLOR® Nichtionische Tenside 15

Test 0-47

Nonionic surfactants / Tensio-actifs
non ioniques / Tensioactivos no iónicos


REF 985 047

620 nm


Method(e) / Método

0471


0.3 - 15.0 mg/l Triton® X-100

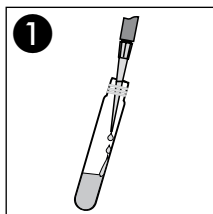


1
4.0 ml Probe
Sample
Enchantillon
Muestra

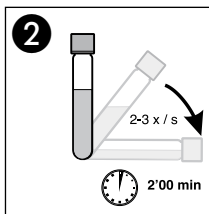


2
Umschwenken
Shake gently
Secouer légèrement
Mezclar volteándolo

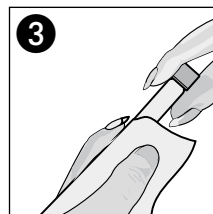

3
Säubern
Clean
Nettoyer
Limpiar

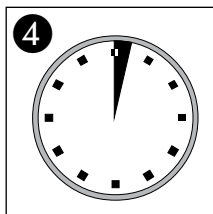

4
2'00 min
Phasentrennung abwarten
wait for phase separation
laisser séparer les phases
esperar la separación de fases

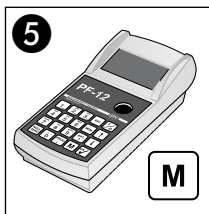

5
Messung
Measurement
Mesure
Medición



Meerwasser
Sea water
Eau de mer
Agua de mar