

NANOCOLOR® Silber 3

Silver / Argent / Plata

620 nm


Method(e) / Método

0491


0.20 - 3.00 mg/l Ag⁺

Test 0-49


REF 985 049

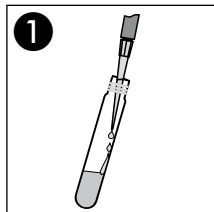


500 µl R2

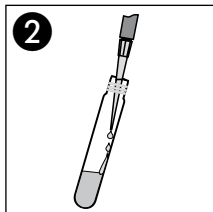


4.0 ml Probe
Sample
Enchantillon
Muestra

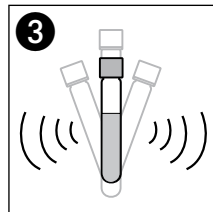


Schütteln
Shake
Agiter
Agitar

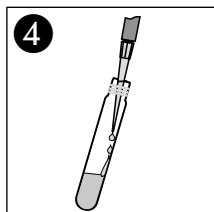


500 µl R3

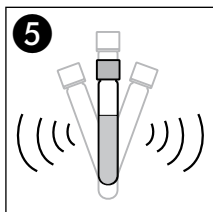


Schütteln
Shake
Agiter
Agitar

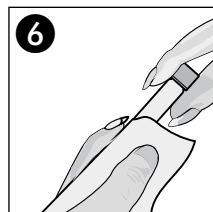

Säubern
Clean
Nettoyer
Limpiar

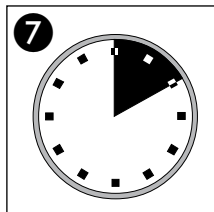

10'00 min

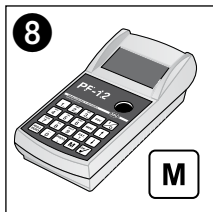

Messung
Measurement
Mesure
Medición

Meerwasser
Sea water
Eau de mer
Agua de mar

