

NANOCOLOR® organische Säuren 3000

organic Acids / Acides organiques
Acidos orgánicos


Test 0-50

REF 985 050


470 nm

Method(e) / Método


0501	30 - 3000	mg/l CH ₃ COOH
0502	0.5 - 50.0	mmol/l CH ₃ COOH

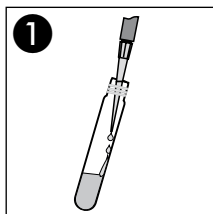


500 µl R2

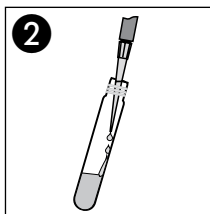


1.0 ml Probe
Sample
Enchantillon
Muestra

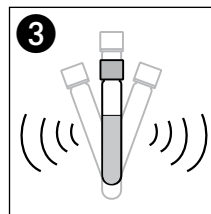


Schütteln
Shake
Agiter
Agitar

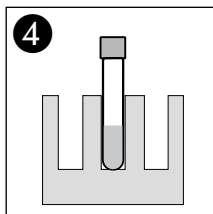


100 °C / 10 min

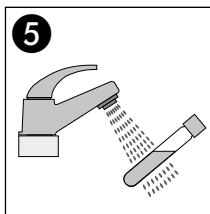


Sofort abkühlen
Cool **immediately**
Refroidir **immédiatement**
Enfriar **inmediata**

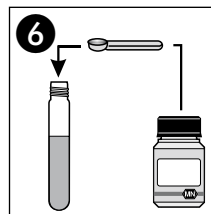


1 x R3

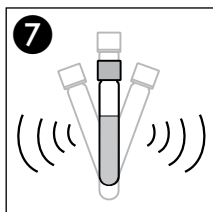


Schütteln
Shake
Agiter
Agitar

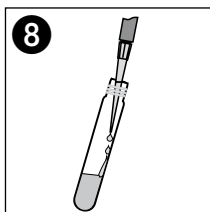


1.0 ml R4

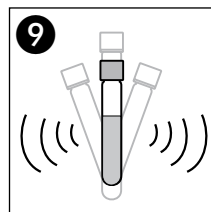


Schütteln
Shake
Agiter
Agitar

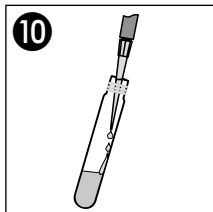


2.0 ml R5

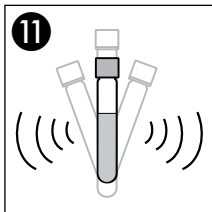


Schütteln
Shake
Agiter
Agitar

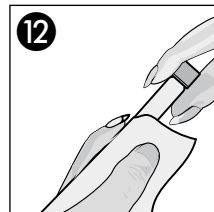


Säubern
Clean
Nettoyer
Limpiar

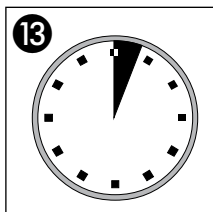

3'00 min

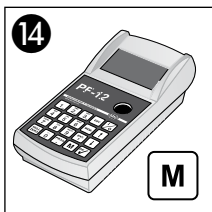

Messung
Measurement
Mesure
Medición

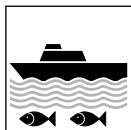

Meerwasser
Sea water
Eau de mer
Agua de mar