

NANOCOLOR® Org. Komplexbildner 10

Complexing agents / Complexants organiques
 Agente formador de complejos org.

Test 0-52

REF 985 052

540 nm


Method(e) / Método

0521


0.5 - 10.0 mg/l I_{Bik}

Nullwert / Blanc value / Zéro / Cero


Messwert / Sample / Echantillon / Muestra

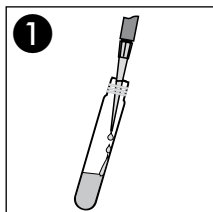


4.0 ml dest. Wasser
 dist. water
 d'eau distillée
 agua dest.

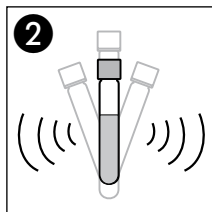


Schütteln
 Shake
 Agiter
 Agitar

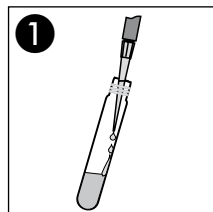


4.0 ml Probe
 Sample
 Echantillon
 Muestra

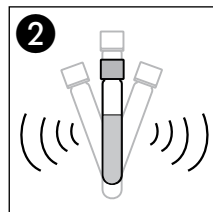


Schütteln
 Shake
 Agiter
 Agitar

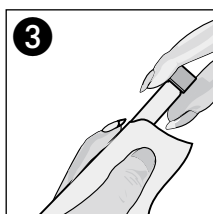

Säubern
 Clean
 Nettoyer
 Limpiar

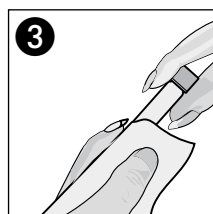


Säubern
 Clean
 Nettoyer
 Limpiar

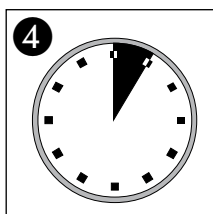


5'00 min




Messung
 Measurement
 Mesure
 Medición

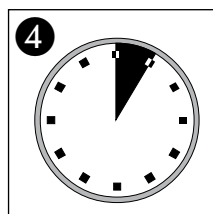

5'00 min

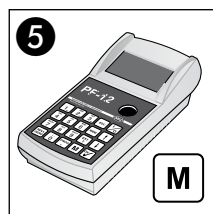

Messung
 Measurement
 Mesure
 Medición

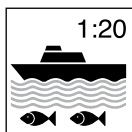

Meerwasser / Sea water
 Eau de mer / Agua de mar