

NANOCOLOR® Kupfer 7

Copper / Cuivre / Cobre

585 nm


Method(e) / Método

0541


0.10 - 7.00 mg/l Cu²⁺

Test 0-54


REF 985 054

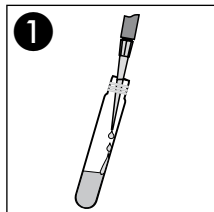


1
4.0 ml Probe
Sample
Enchantillon
Muestra

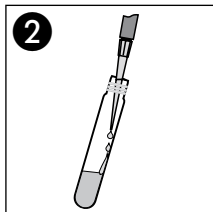


2
200 µl R2

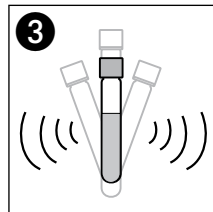


3
Schütteln
Shake
Agiter
Agitar

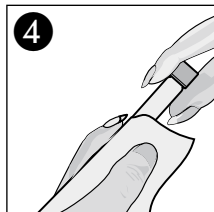


4
Säubern
Clean
Nettoyer
Limpiar

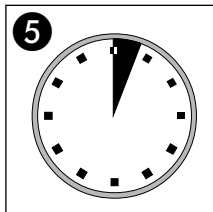

5
10'00 min

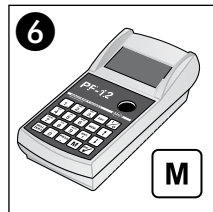

6
Messung
Measurement
Mesure
Medición

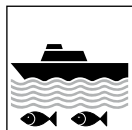

Meerwasser
Sea water
Eau de mer
Agua de mar