

# NANOCOLOR® Mangan 10

Manganese / Manganèse / Manganeso

**Test 0-58**


REF 985 058

470 nm


Method(e) / Método

0581


0.1 - 10.0 mg/l Mn

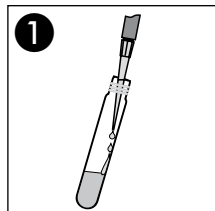


**1**  
4.0 ml Probe  
Sample  
Enchantillon  
Muestra

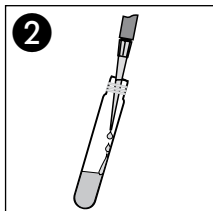


**2**  
500 µl R2

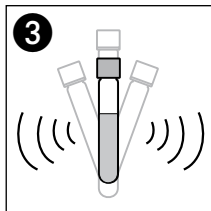


**3**  
Schütteln  
Shake  
Agiter  
Agitar

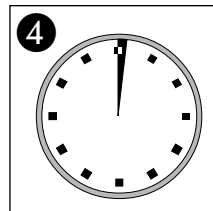


**4**  
1'00 min

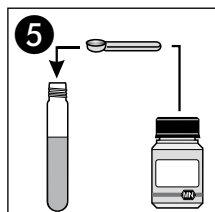


**5**  
1 x R3

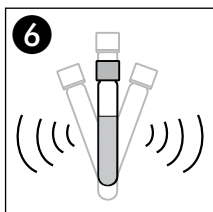


**6**  
Kräftig Schütteln  
Shake well  
Bien Agiter  
Agitar intensamente

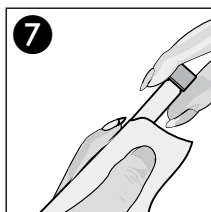


**7**  
Säubern  
Clean  
Nettoyer  
Limpiar

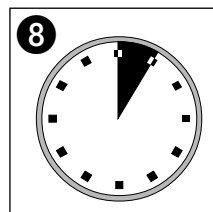

**8**  
5'00 min

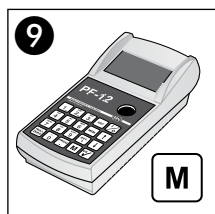

Messung  
Measurement  
Mesure  
Medición

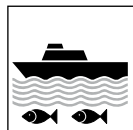

Meerwasser  
Sea water  
Eau de mer  
Agua de mar