

NANOCOLOR® Nickel 7

Níquel

470 nm


Method(e) / Método

0611


0.10 - 7.00 mg/l Ni²⁺

Test 0-61


REF 985 061

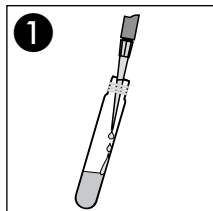


1
5.0 ml Probe
Sample
Enchantillon
Muestra

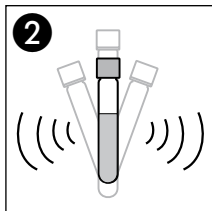


2
Schütteln
Shake
Agiter
Agitar

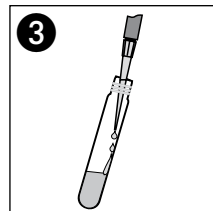


3
1.0 ml R2

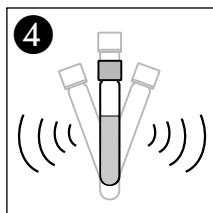


4
Schütteln
Shake
Agiter
Agitar

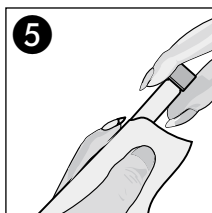


5
Säubern
Clean
Nettoyer
Limpiar

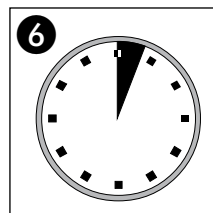

6
3'00 min

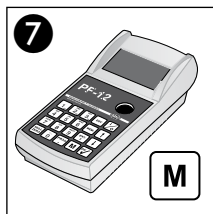

7
Messung
Measurement
Mesure
Medición

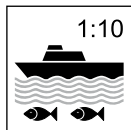

1:10

Meerwasser
Sea water
Eau de mer
Agua de mar