

NANOCOLOR® Nitrat 50

Nitrate / Nitrato


345 nm

Method(e) / Método


0641 0.3 - 22.0 mg/l NO₃-N
0642 2 - 100 mg/l NO₃⁻

Test 0-64


REF 985 064

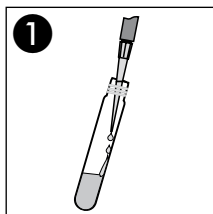


500 µl Probe
Sample
Enchantillon
Muestra

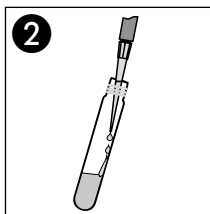


500 µl R2

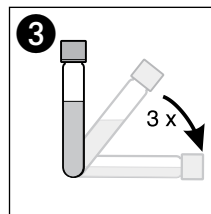


Umschwenken
Shake gently
Secouer légèrement
Mezclar volteándolo

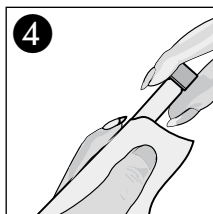

Säubern
Clean
Nettoyer
Limpiar

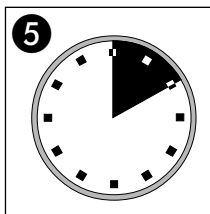

10'00 min

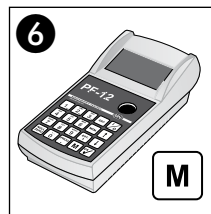

Messung
Measurement
Mesure
Medición



Meerwasser
Sea water
Eau de mer
Agua de mar