

NANOCOLOR® Nitrat 8

Nitrate / Nitrato


345 nm

Method(e) / Método


0651 0.30 - 8.00 mg/l NO₃-N
0652 1.3 - 35.0 mg/l NO₃⁻

Test 0-65


REF 985 065

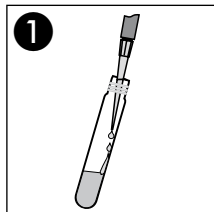


1
500 µl Probe
Sample
Echantillon
Muestra

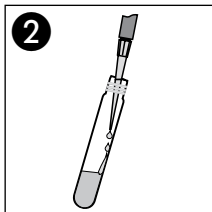


2
500 µl R2

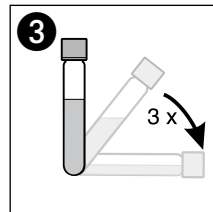


3
Umschwenken
Shake gently
Secouer légèrement
Mezclar volteándolo

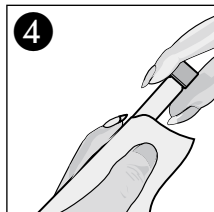

4
Säubern
Clean
Nettoyer
Limpiar

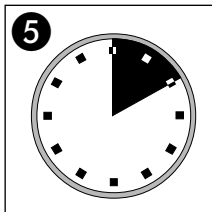

5
10'00 min

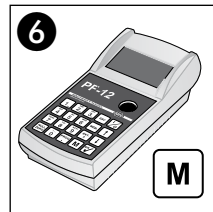

6
Messung
Measurement
Mesure
Medición



Meerwasser
Sea water
Eau de mer
Agua de mar