

NANOCOLOR® Nitrit 2

Nitrite / Nitrito


540 nm

Method(e) / Método


0681 0.003 - 0.460 mg/l NO₂-N
0682 0.02 - 1.50 mg/l NO₂⁻

Test 0-68


REF 985 068

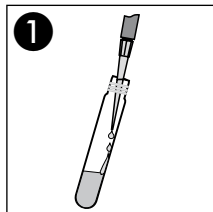


1
4.0 ml Probe
Sample
Enchantillon
Muestra

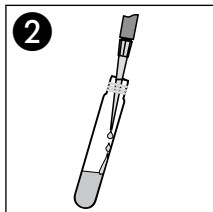


2
500 µl R2

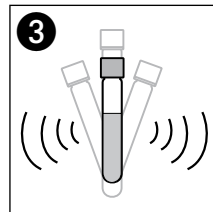


3
Schütteln
Shake
Agiter
Agitar

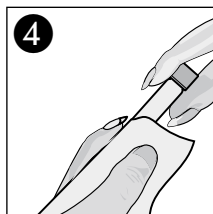


4
Säubern
Clean
Nettoyer
Limpiar

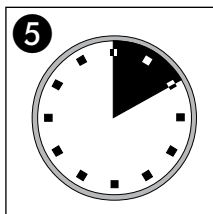

5
10'00 min

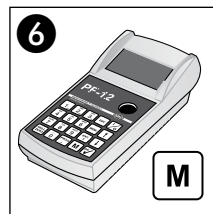

6
Messung
Measurement
Mesure
Medición

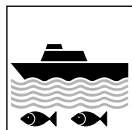

Meerwasser
Sea water
Eau de mer
Agua de mar