

NANOCOLOR® Sulfid 3

Sulphide / Sulfure / Sulfuro

620 nm


Method(e) / Método

0731


0.05 - 3.00 mg/l S²⁻

Test 0-73


REF 985 073

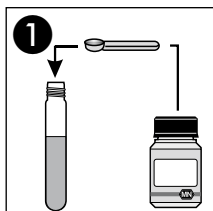


1 x R2

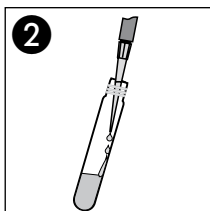


4.0 ml Probe
Sample
Echantillon
Muestra

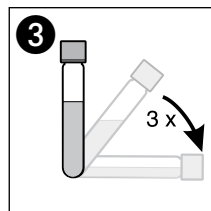


Umschwenken
Shake gently
Secouer légèrement
Mezclar volteándolo

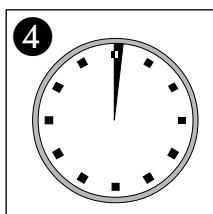


1'00 min

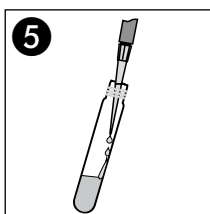


200 µl R3

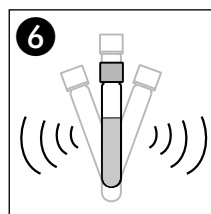


Schütteln
Shake
Agiter
Agitar

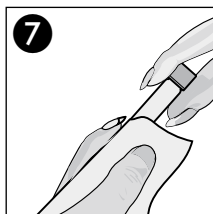


Säubern
Clean
Nettoyer
Limpiar

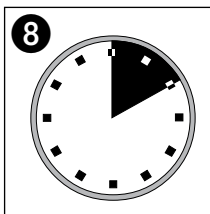

10'00 min

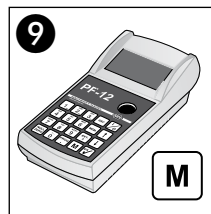

Messung
Measurement
Mesure
Medición

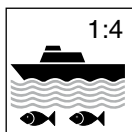

1:4
Meerwasser
Sea water
Eau de mer
Agua de mar