

# NANOCOLOR® Resthärte 1

Residual Hardness / Dureté résiduelle  
Dureza residual

540 nm


Method(e) / Método

0841/2 0.02 - 1.00 °d / 0.004 - 0.180 mmol/l


0843/4 0.03 - 1.25 °e / 0.04 - 1.78 °f

**Test 0-84**


REF 985 084

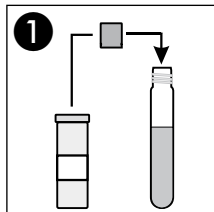


1 x NANOFIX R2

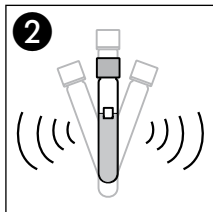


Kräftig Schütteln  
Shake well  
Bien Agiter  
Agitar intencamente

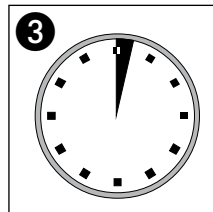


2'00 min

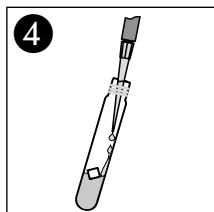


5.0 ml Probe  
Sample  
Enchantillon  
Muestra

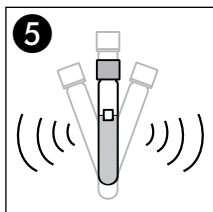


Schütteln  
Shake  
Agiter  
Agitar

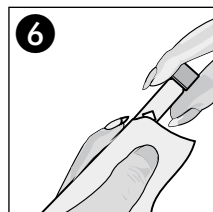

Säubern  
Clean  
Nettoyer  
Limpiar

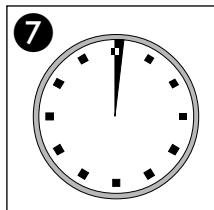

1'00 min

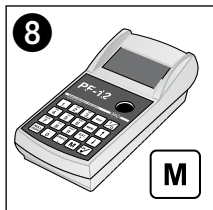

Messung  
Measurement  
Mesure  
Medición

Meerwasser  
Sea water  
Eau de mer  
Agua de mar

