

NANOCOLOR® Stärke 100

Starch / Amidon / Almidó

540 nm


Method(e) / Método

0851


5 - 100 mg/l

Test 0-85


REF 985 085

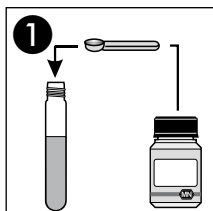


1 x R2

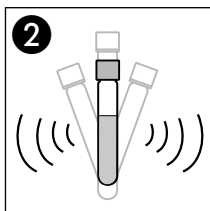


Schütteln
Shake
Agiter
Agitar

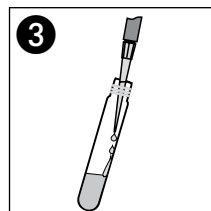


1.0 ml R3

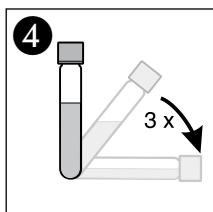


Umschwenken
Shake gently
Secouer légèrement
Mezclar volteándolo

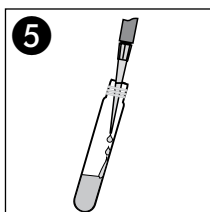


4.0 ml Probe
Sample
Enchantillon
Muestra

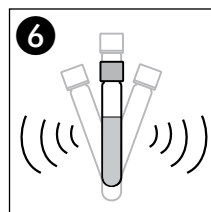


Schütteln
Shake
Agiter
Agitar

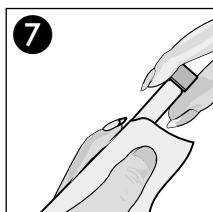


Säubern
Clean
Nettoyer
Limpiar

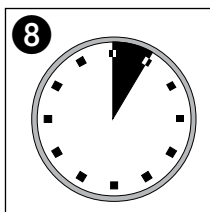

5'00 min

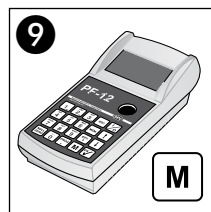

Messung
Measurement
Mesure
Medición

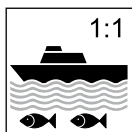

Meerwasser
Sea water
Eau de mer
Agua de mar