

# NANOCOLOR® Sulfat 200

Sulphate / Sulfate / Sulfato

436 nm


Method(e) / Método

0861


10 - 200 mg/l  $\text{SO}_4^{2-}$

**Test 0-86**


REF 985 086

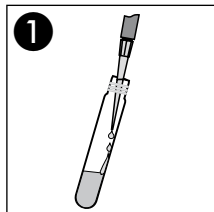


**1**  
4.0 ml Probe  
Sample  
Enchantillon  
Muestra

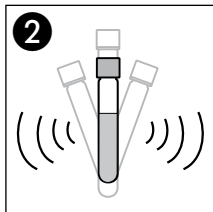


**2**  
Schütteln  
Shake  
Agiter  
Agitar

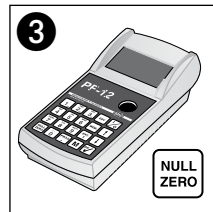


**3**  
Nullwert  
Blank value  
Zéro  
Cero

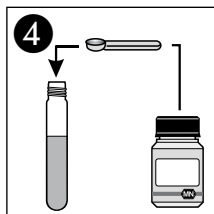


**4**  
1 x R2

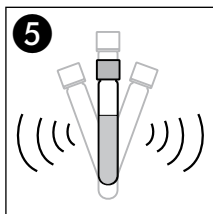


**5**  
10 s Schütteln  
Shake  
Agiter  
Agitar

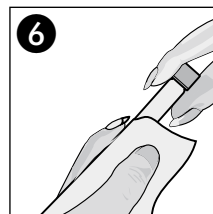

**6**  
Säubern  
Clean  
Nettoyer  
Limpiar

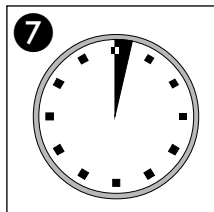

**7**  
2'00 min

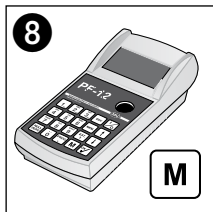

**8**  
Messung  
Measurement  
Mesure  
Medición

Meerwasser  
Sea water  
Eau de mer  
Agua de mar

