

NANOCOLOR® Sulfat 1000

Sulphate / Sulfate / Sulfato

436 nm


Method(e) / Método

0871


200 - 1000 mg/l SO_4^{2-}

Test 0-87

REF 985 087

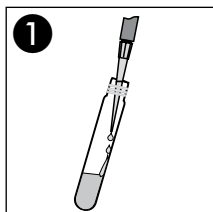


1
1.0 ml Probe
Sample
Enchantillon
Muestra

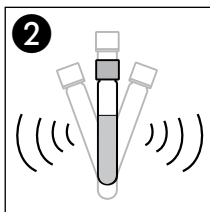


2
Schütteln
Shake
Agiter
Agitar




3
Nullwert
Blank value
Zéro
Cero

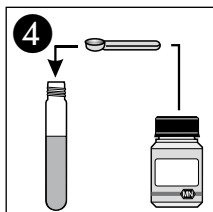


4
1 x R2

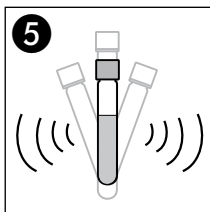


5
10 s Schütteln
Shake
Agiter
Agitar

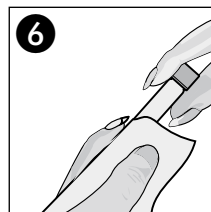

6
Säubern
Clean
Nettoyer
Limpiar

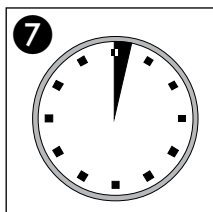

7
2'00 min

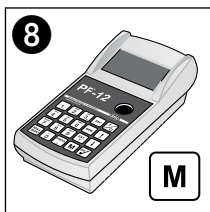

8
Messung
Measurement
Mesure
Medición

Meerwasser
Sea water
Eau de mer
Agua de mar

