

# NANOCOLOR® Sulfit 10

Sulphite / Sulfite / Sulfito

436 nm


Method(e) / Método

0891


0.2 - 10.0 mg/l  $\text{SO}_3^{2-}$

**Test 0-89**


REF 985 089

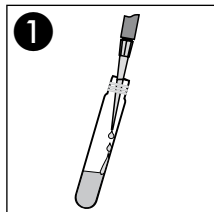


**1**  
4.0 ml Probe  
Sample  
Enchantillon  
Muestra

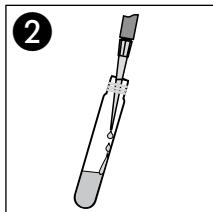


**2**  
200 ml R2

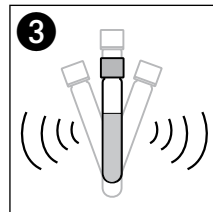


**3**  
Schütteln  
Shake  
Agiter  
Agitar

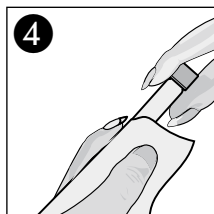


**4**  
Säubern  
Clean  
Nettoyer  
Limpiar

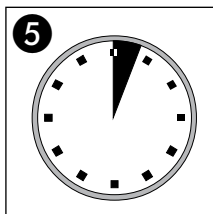

**5**  
5'00 min

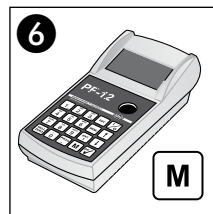

**6**  
Messung  
Measurement  
Mesure  
Medición

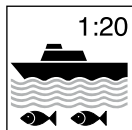

1:20

Meerwasser  
Sea water  
Eau de mer  
Agua de mar