

# NANOCOLOR® Sulfit 100

Sulphite / Sulfite / Sulfito

470 nm

Method(e) / Método


0901

5 - 100 mg/l  $\text{SO}_3^{2-}$


**Test 0-90**

REF 985 090


## Nullwert / Blanc value / Zéro / Cero

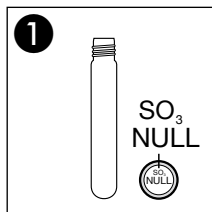


**Leere Rundküvette**  
**Empty test tube**  
Cuve ronde **vide**  
Tubo de test **vacío**

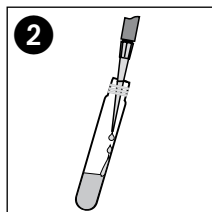

1.0 ml R1

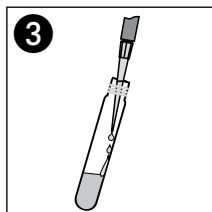


200 µl R2

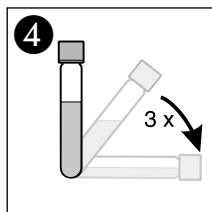


Umschwenken  
Shake gently  
Secouer légèrement  
Mezclar volteándolo


## Messwert / Sample / Echantillon / Muestra

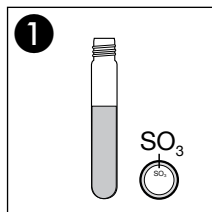

**Sulfit-Rundküvette**  
Test tube **Sulphite**  
Cuve ronde **Sulfite**  
Tubo de test **Sulfito**

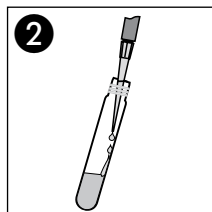


200 µl R2

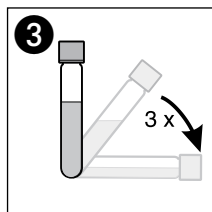


Umschwenken  
Shake gently  
Secouer légèrement  
Mezclar volteándolo


**Nullwert / Blanc value / Zéro / Cero**

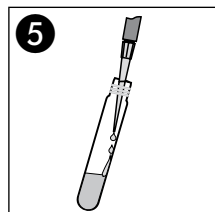


**4.0 ml Probe**  
Sample  
Echantillon  
Muestra

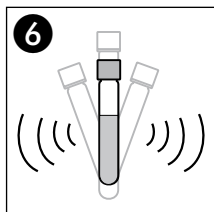


**Schütteln**  
Shake  
Agiter  
Agitar

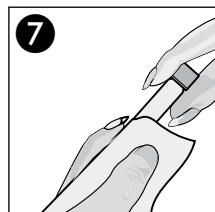

**Säubern**  
Clean  
Nettoyer  
Limpiar

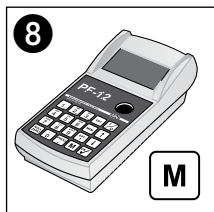


**Messung**  
Measurement  
Mesure  
Medición

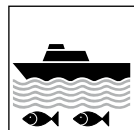


**Meerwasser**  
Sea water  
Eau de mer  
Agua de mar


**Messwert / Sample / Echantillon / Muestra**

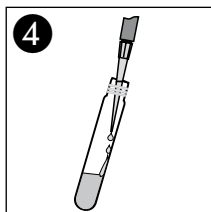


**4.0 ml Probe**  
Sample  
Echantillon  
Muestra

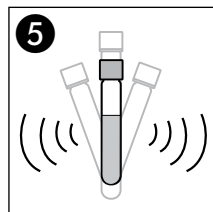

**Schütteln**  
Shake  
Agiter  
Agitar

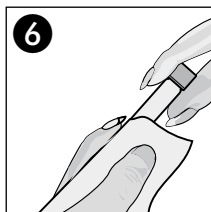

**Säubern**  
Clean  
Nettoyer  
Limpiar

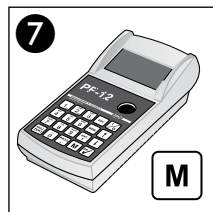

**Messung**  
Measurement  
Mesure  
Medición