

NANOCOLOR® TOC 25 COT

585 nm


Method(e) / Método

0931


2.0 - 25.0 mg/l C

Test 0-93


REF 985 093

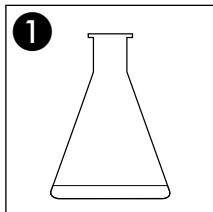


1
10.0 ml Probe
Sample
Enchantillon
Muestra

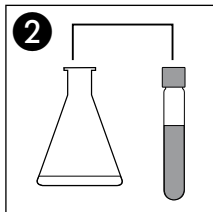


2
500 µl R0

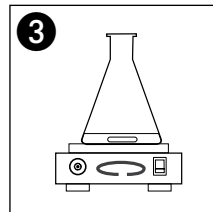


3
10 min Rühren
Stirr
Agiter
Agitar

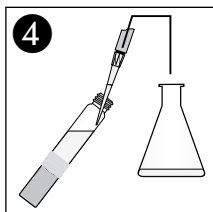


4
5.0 ml in TOC-Küvette
into tube TOC
dans cuve COT
en tubo TOC

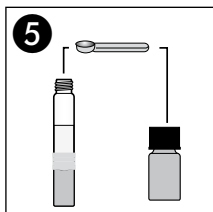


5
1 x R2

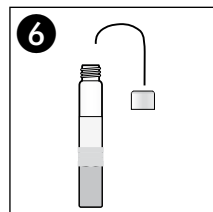


6
Thermokappe aufschrauben
Close with thermo cap
Fermer avec thermo cap
Cerrar con thermo cap

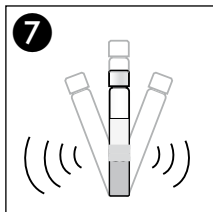


7
Schütteln
Shake
Agiter
Agitar

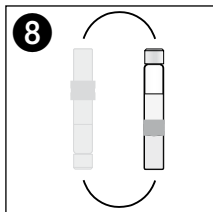


8
Küvette um 180° drehen
Turn tube upside down
Retourner tube
Girar tubo

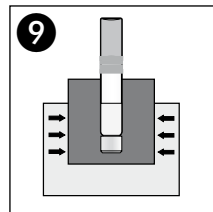


9
120 °C / 2 h

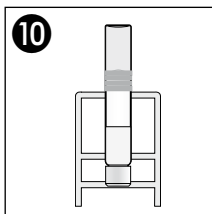


60 min Abkühlen lassen
Allow to cool
Laisser refroidir
Dejar enfriar

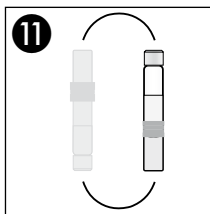

Küvette um 180° drehen
Turn tube upside down
Retourner tube
Girar tubo

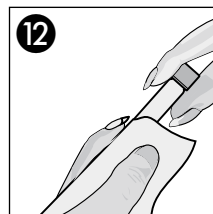

Säubern
Clean
Nettoyer
Limpiar

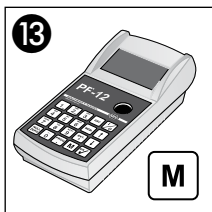

Messung
Measurement
Mesure
Medición



Meerwasser
Sea water
Eau de mer
Agua de mar