

# NANOCOLOR® Aluminium 07

Aluminio

540 nm

Method(e) / Método

0981


0.02 - 0.70 mg/l Al<sup>3+</sup>

**Test 0-98**


REF 985 098




Null Messung  
Measure blank  
Lecture blanc  
Lectura blanco

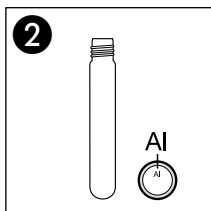


Aluminium-Rundküvette  
Test tube Aluminium  
Cuve ronde Aluminium  
Tubo de test Aluminio

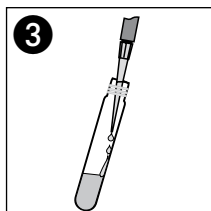


500 µl R2

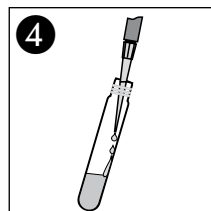


500 µl R3

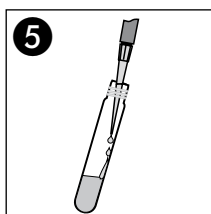


4.0 ml Probe  
Sample  
Enchantillon  
Muestra

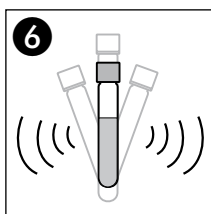


Schütteln  
Shake  
Agiter  
Agitar

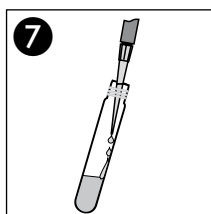


500 µl R4

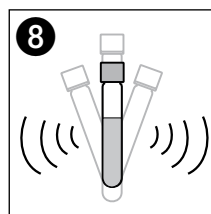


Schütteln  
Shake  
Agiter  
Agitar

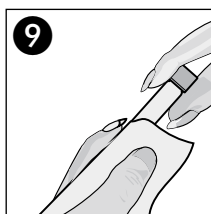


Säubern  
Clean  
Nettoyer  
Limpiar

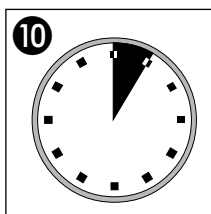

5'00 min

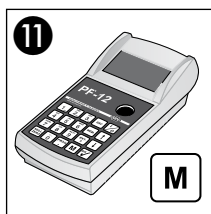

Messung  
Measurement  
Mesure  
Medición

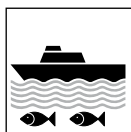

Meerwasser  
Sea water  
Eau de mer  
Agua de mar