

NANOCOLOR® Methanol 15

Méthanol / Metanol

Test 8-59


REF 985 859

620 nm


Method(e) / Método

0591


0.2 - 15.0 mg/l MeOH

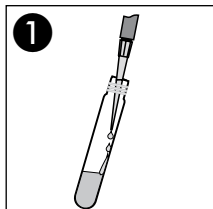


3.0 ml R1

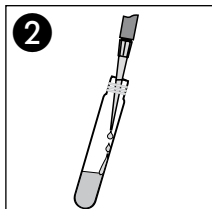


1.5 ml Probe
Sample
Enchantillon
Muestra

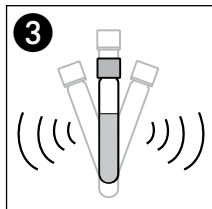


Schütteln
Shake
Agiter
Agitar

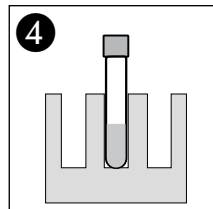


25 °C / 30 min

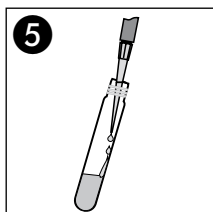


100 µl R2

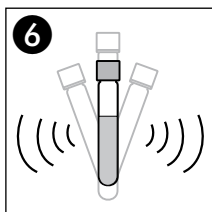


Schütteln
Shake
Agiter
Agitar

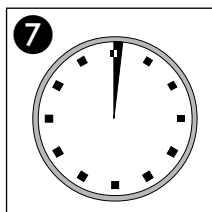


1'00 min

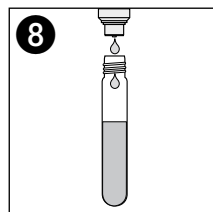


2 x R3

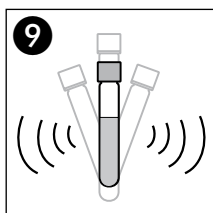


Schütteln
Shake
Agiter
Agitar

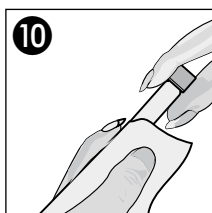

Säubern
Clean
Nettoyer
Limpiar

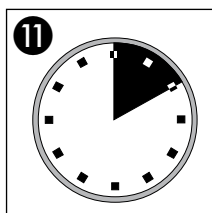

10'00 min

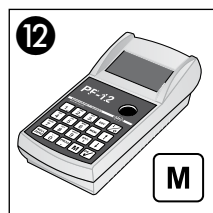

Messung
Measurement
Mesure
Medición



Meerwasser
Sea water
Eau de mer
Agua de mar