

visocolor[®] Eisen (Triazin) Iron (Triazine) / Fer (Triazine) / Hierro (Triacina)

540 nm


Method(e) / Método

5251


0.10 - 2.00 mg/l Fe

Test 5-25


REF 914 239

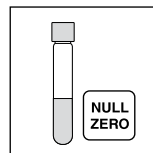


Probe / sample
échantillon / muestra

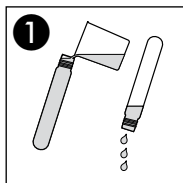


1
Küvetten spülen
Rinse test tubes
Rincer les éprouvettes
Lavar los tubos

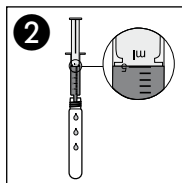


2
5 ml Probe
Sample
Echantillon
Muestra

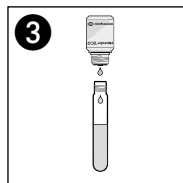


3
5 Δ Eisen-1
Iron-1
Fer-1
Hierro-1

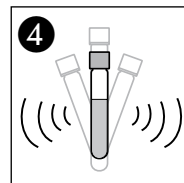


4
Schütteln
Shake
Agiter
Agitar

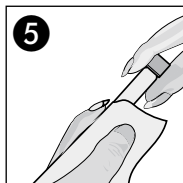


5
Säubern
Clean
Nettoyer
Limpiar

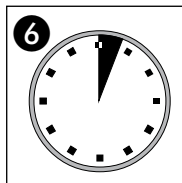

6
3'00 min

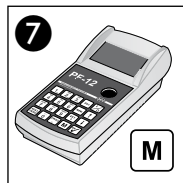

7
Messung
Measurement
Mesure
Medición

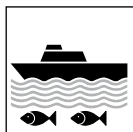

Meerwasser
Sea water
Eau de mer
Agua de mar