

# visicolor® ECO Brom (frei)

Bromine (free) / Brome (libre) / Bromo (libre)

540 nm


Method(e) / Método

5111


0.10–13.00 mg/L Br<sub>2</sub>

**Test 5-11**


REF 931 211

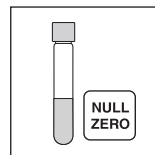


Probe / sample  
échantillon / muestra

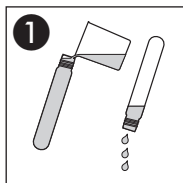


**1**  
Küvetten spülen  
Rinse test tubes  
Rincer les éprouvettes  
Lavar los tubos

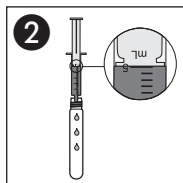


**2**  
5 mL Probe  
Sample  
Echantillon  
Muestra

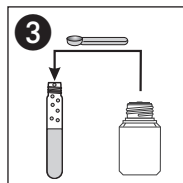


**3** 1 Br<sub>2</sub>-1

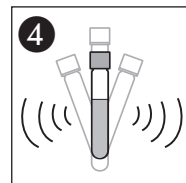


**4**  
20 s Schütteln  
Shake  
Agiter  
Agitar

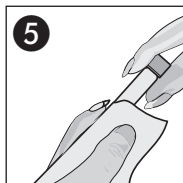

**5**  
Säubern  
Clean  
Nettoyer  
Limpiar

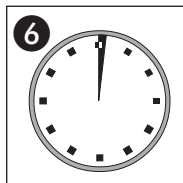

**6**  
1'00 min

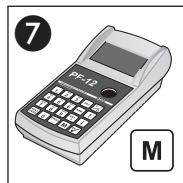

**7**  
Messung  
Measurement  
Mesure  
Medición



Meerwasser  
Sea water  
Eau de mer  
Agua de mar

# visicolor® ECO Brom (gesamt)

Bromine (total) / Brome (totale) / Bromo (total)

540 nm


Method(e) / Método

5112


0.10–13.00 mg/L Br<sub>2</sub>

**Test 5-11**


REF 931 211

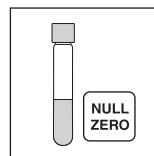


Probe / sample  
échantillon / muestra

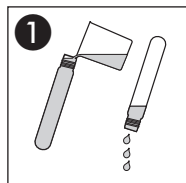


**1**  
Küvetten spülen  
Rinse test tubes  
Rincer les éprouvettes  
Lavar los tubos

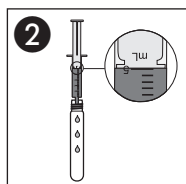


**2**  
5 mL Probe  
Sample  
Echantillon  
Muestra

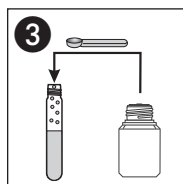


**3**  
1 Pipette Br<sub>2</sub>-1

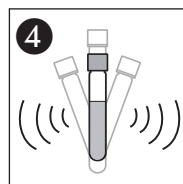


**4**  
20 s Schütteln  
Shake  
Agiter  
Agitar

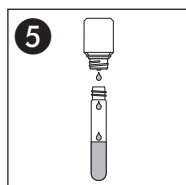


**5**  
3 Drops Br<sub>2</sub>-2

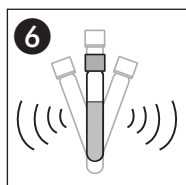


**6**  
Schütteln  
Shake  
Agiter  
Agitar

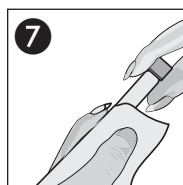

**7**  
Säubern  
Clean  
Nettoyer  
Limpiar

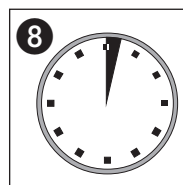

**8**  
2'00 min

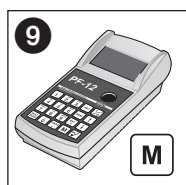

Messung  
Measurement  
Mesure  
Medición



Meerwasser  
Sea water  
Eau de mer  
Agua de mar