

visicolor® ECO Chlor 6 (frei)

Chlorine (free) / Chlore (libre) / Cloro (libre)


Test 5-17/5-19

REF 931 217/931 219


540 nm

Method(e) / Método


5171/5191 0.05–6.00 mg/L Cl₂

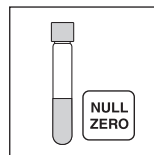


Probe / sample
échantillon / muestra

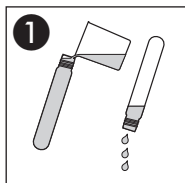


1
Küvetten spülen
Rinse test tubes
Rincer les éprouvettes
Lavar los tubos

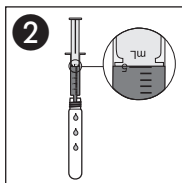


2
5 mL Probe
Sample
Echantillon
Muestra

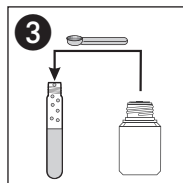


3
1  Cl₂-1

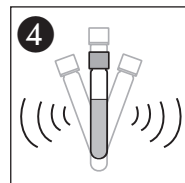


4
20 s Schütteln
Shake
Agiter
Agitar

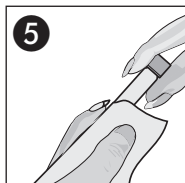

5
Säubern
Clean
Nettoyer
Limpiar

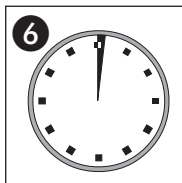

6
1'00 min

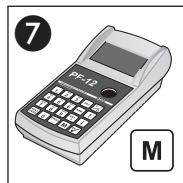

7
Messung
Measurement
Mesure
Medición



Meerwasser
Sea water
Eau de mer
Agua de mar