

visicolor® ECO Eisen 2

Iron / Fer / Hierro

540 nm


Method(e) / Método

5261


0.04–2.00 mg/L Fe

Test 5-26


REF 931 226

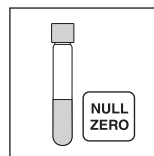


Probe / sample
échantillon / muestra

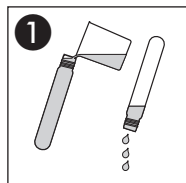


1
Küvetten spülen
Rinse test tubes
Rincer les éprouvettes
Lavar los tubos

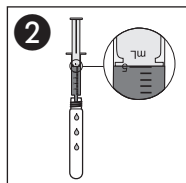


2
5 mL Probe
Sample
Echantillon
Muestra

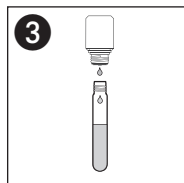


3
4 ♀ Fe-1

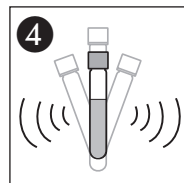


4
Schütteln
Shake
Agiter
Agitar

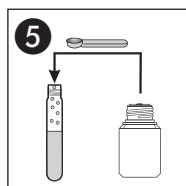


5
1 Fe-2

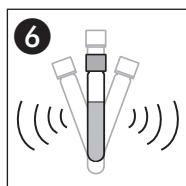


6
Schütteln
Shake
Agiter
Agitar

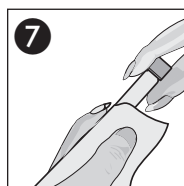

7
Säubern
Clean
Nettoyer
Limpiar

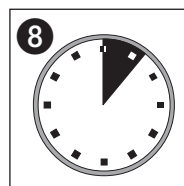

8
7'00 min

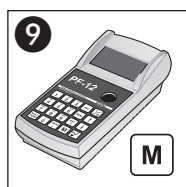

Messung
Measurement
Mesure
Medición



Meerwasser
Sea water
Eau de mer
Agua de mar