

visicolor® ECO Kupfer

Copper / Cuivre / Cobre

585 nm


Method(e) / Método

5371


0.1–5.0 mg/L Cu²⁺

Test 5-37


REF 931 237

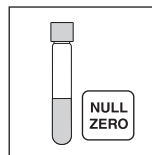


Probe / sample
échantillon / muestra

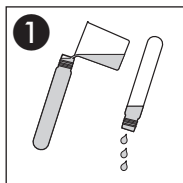


1
Küvetten spülen
Rinse test tubes
Rincer les éprouvettes
Lavar los tubos

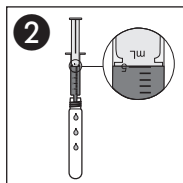


2
5 mL Probe
Sample
Echantillon
Muestra

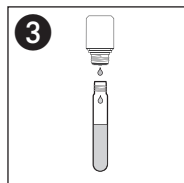


3
5 Δ Cu-1

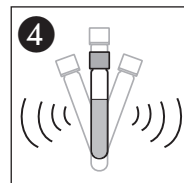


4
Schütteln
Shake
Agiter
Agitar

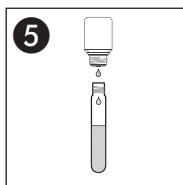


5
5 Δ Cu-2

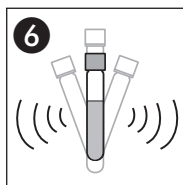


6
Schütteln
Shake
Agiter
Agitar

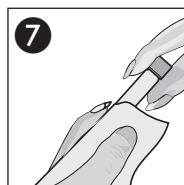

7
Säubern
Clean
Nettoyer
Limpiar

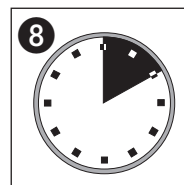

8
10'00 min

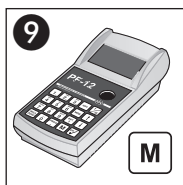

Messung
Measurement
Mesure
Medición



Meerwasser
Sea water
Eau de mer
Agua de mar