

visicolor® ECO Mangan

Manganese / Manganèse / Manganeso

436 nm


Method(e) / Método

5381


0.1–5.0 mg/L Mn

Test 5-38


REF 931 238

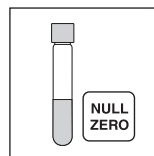


Probe / sample
échantillon / muestra

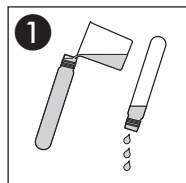


1
Küvetten spülen
Rinse test tubes
Rincer les éprouvettes
Lavar los tubos

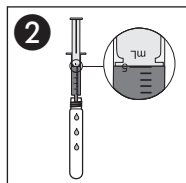


2
5 mL Probe
Sample
Echantillon
Muestra

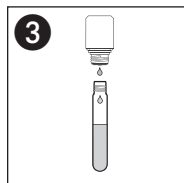


3
5 ∆ Mn-1

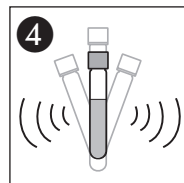


4
Schütteln
Shake
Agiter
Agitar

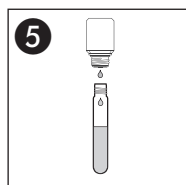


5
7 ∆ Mn-2

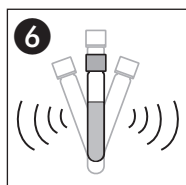


6
Schütteln
Shake
Agiter
Agitar

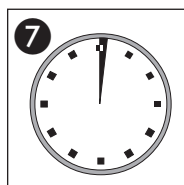


7
1'00 min

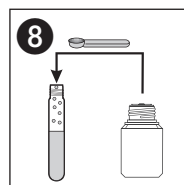


8
1 ∆ Mn-3

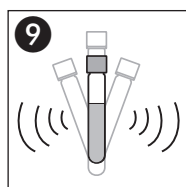


9
Schütteln
Shake
Agiter
Agitar

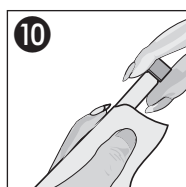

10
Säubern
Clean
Nettoyer
Limpiar

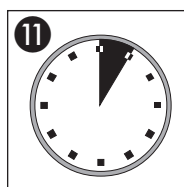

11
5'00 min

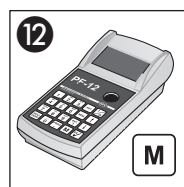

12
Messung
Measurement
Mesure
Medición



Meerwasser
Sea water
Eau de mer
Agua de mar