

visicolor® ECO Nitrit

Nitrite / Nitrito

540 nm

Method(e) / Método

5441


0.01–0.15 mg/L NO₂⁻N

5442


0.02–0.50 mg/L NO₂⁻

Test 5-44


REF 931 244

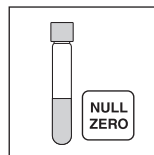


Probe / sample
échantillon / muestra

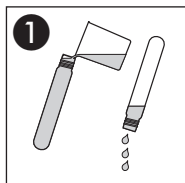


Küvetten spülen
Rinse test tubes
Rincer les éprouvettes
Lavar los tubos

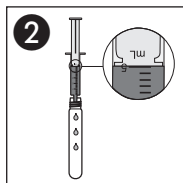


5 mL Probe
Sample
Echantillon
Muestra

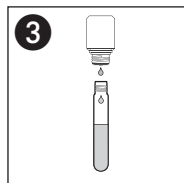


4 Δ NO₂-1

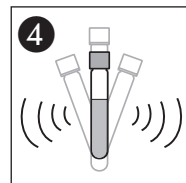


Schütteln
Shake
Agiter
Agitar

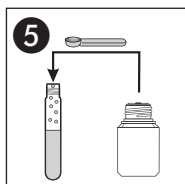


1 Δ NO₂-2

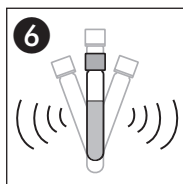


Schütteln
Shake
Agiter
Agitar

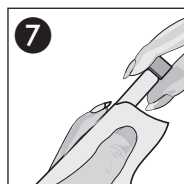

Säubern
Clean
Nettoyer
Limpiar

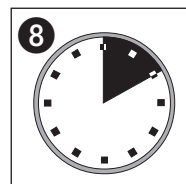

10'00 min

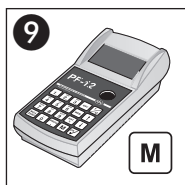

Messung
Measurement
Mesure
Medición



Meerwasser
Sea water
Eau de mer
Agua de mar