

visicolor® ECO pH 6.0–8.2

Test 5-70


REF 931 270

436/540 nm


Method(e) / Método

5701


6.0–8.2 pH

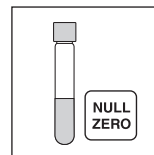


Probe / sample
échantillon / muestra

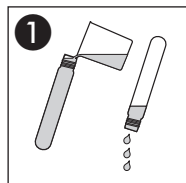


1
Küvetten spülen
Rinse test tubes
Rincer les éprouvettes
Lavar los tubos

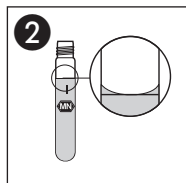


2
10 mL Probe
Sample
Echantillon
Muestra

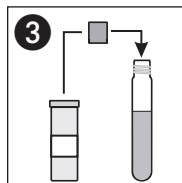


3
1 NANOPIX pH

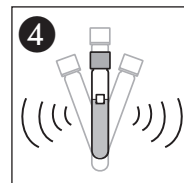


4
Schütteln
Shake
Agiter
Agitar

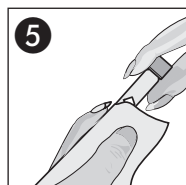

5
Säubern
Clean
Nettoyer
Limpiar

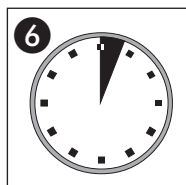

6
3'00 min

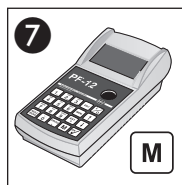

7
Messung
Measurement
Mesure
Medición



Meerwasser
Sea water
Eau de mer
Agua de mar