

# visicolor® ECO Phosphat

Phosphate / Fosfato

**Test 5-84**

REF 931 284

690 nm


Method(e) / Método

5841


0.2–5.0 mg/L PO<sub>4</sub>-P

5842


0.6–15.0 mg/L PO<sub>4</sub><sup>3-</sup>

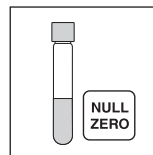


Probe / sample  
échantillon / muestra

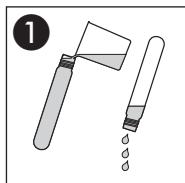


Küvetten spülen  
Rinse test tubes  
Rincer les éprouvettes  
Lavar los tubos

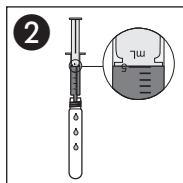


5 mL Probe  
Sample  
Echantillon  
Muestra

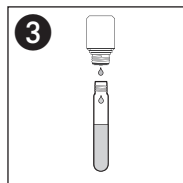


6 Δ PO<sub>4</sub>-1

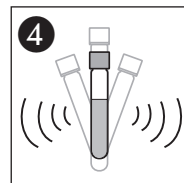


Schütteln  
Shake  
Agiter  
Agitar

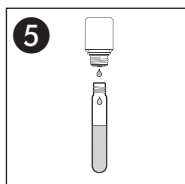


6 Δ PO<sub>4</sub>-2

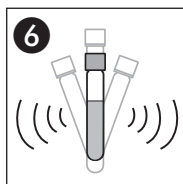


Schütteln  
Shake  
Agiter  
Agitar

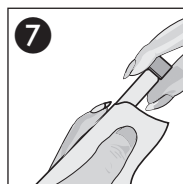

Säubern  
Clean  
Nettoyer  
Limpiar

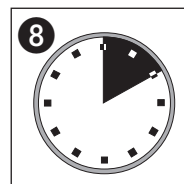

10'00 min

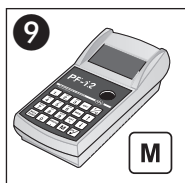

Messung  
Measurement  
Mesure  
Medición



Meerwasser  
Sea water  
Eau de mer  
Agua de mar