

# Zirconium 100

Zirconium 100

Circonio 100

Rundküvettestest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica

5–100 mg/L Zr

**Test 0-01**


REF 985 001




540 nm

Method(e) / Método


0011

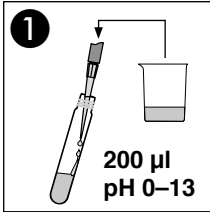


4.0 mL Probe  
Sample  
Echantillon  
Muestra

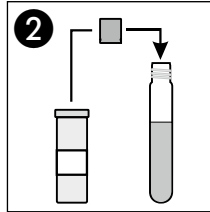


1 x **NANOFIX R2**

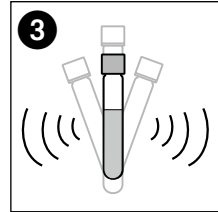


Schütteln  
Shake  
Agiter  
Agitar

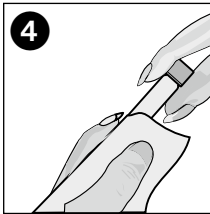

Säubern  
Clean  
Nettoyer  
Limpiar

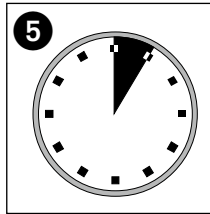

5'00 min

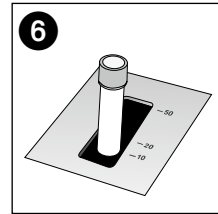

Messung  
Measurement  
Mesure  
Medición