

Ammonium 2000

Amonio

Test 0-02

REF 985 002



Rundküvettentest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica

300–1600 mg/L $\text{NH}_4\text{-N}$

400–2000 mg/L NH_4^+

400–2000 mg/L NH_3


585 nm

Method(e) / Método


0021

0022


0023

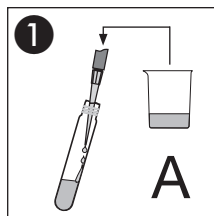


500 μL Probe
Sample
Echantillon
Muestra

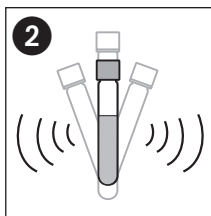


Schütteln
Shake
Agiter
Agitar

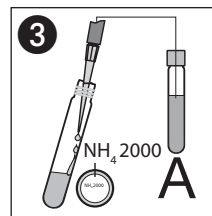


200 μL A

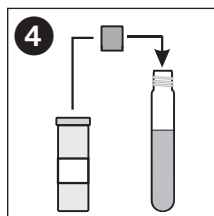


1 x **NANOFIX R2**

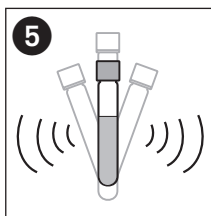


Schütteln
Shake
Agiter
Agitar

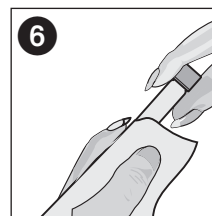

Säubern
Clean
Nettoyer
Limpiar

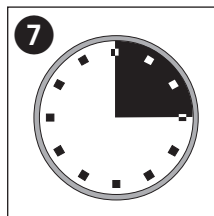

15'00 min

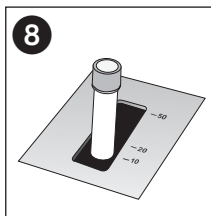

Messung
Measurement
Mesure
Medición